

Evaluation of Lateral-Flow *Clostridium botulinum* Neurotoxin Detection Kits for Food Analysis

Shashi K. Sharma,^{1*} Brian S. Eblen,¹ Robert L. Bull,² Donald H. Burr,³
and Richard C. Whiting¹

U.S. Food and Drug Administration, Center for Food Safety and Applied Nutrition, College Park, Maryland 20740¹;
Biological Defense Research Directorate, Naval Medical Research Center, Silver Spring,
Maryland 20910²; and National Center for Food Safety and Technology,
Summit-Argo, Illinois 60501³

Received 20 April 2004/Accepted 8 February 2005

The suitability and sensitivity of two in vitro lateral-flow assays for detecting *Clostridium botulinum* neurotoxins (BoNTs) in an assortment of foods were evaluated. Toxin extraction and preparation methods for various liquid, solid, and high-fat-content foods were developed. The lateral-flow assays, one developed by the Naval Medical Research Center (Silver Spring, MD) and the other by Alexeter Technologies (Gaithersburg, MD), are based on the immunodetection of BoNT types A, B, and E. The assays were found to be rapid and easy to perform with minimum requirements for laboratory equipment or skills. They can readily detect 10 ng/ml of BoNT types A and B and 20 ng/ml of BoNT type E. Compared to other in vitro detection methods, these assays are less sensitive, and the assessment of a result is strictly qualitative. However, the assay was found to be simple to use and to require minimal training. The assays successfully detected BoNT types A, B, and E in a wide variety of foods, suggesting their potential usefulness as a preliminary screening system for triaging food samples with elevated BoNT levels in the event of a *C. botulinum* contamination event.

Clostridium botulinum produces seven (A to G) structurally related but antigenically distinct protein neurotoxins. These botulinum neurotoxins (BoNTs) are the causative agents of botulism (3, 18, 21, 22, 29). There has been much effort by the food industry to ensure that food treatment processes prevent the growth of *C. botulinum* spores and the formation of their toxin. At present, the mouse bioassay is the only method that can be used with confidence to detect the biologically active BoNTs (9, 11, 12, 14, 20). Although the mouse test is exquisitely sensitive, with a detection limit of 1 mouse lethal dose (MLD), which is equivalent to approximately 10 pg/ml of neurotoxin, it has a number of drawbacks; it is expensive to perform, requires live mice, and is not specific for the neurotoxin serotype unless neutralization tests with a serotype-specific antiserum are carried out in parallel. In addition, the test takes up to 4 days to complete. Because of its disadvantages, there have been numerous attempts to replace the mouse bioassay with an in vitro method.

A number of in vitro immunoassays with sensitivities comparable to that of the mouse bioassay have been described for BoNTs (13, 14, 15, 27, 28). An enzyme-linked coagulation assay was reported to have a sensitivity comparable to that of the mouse bioassay (13). In this assay, a complex consisting of toxin bound to chicken antibody (or biotinylated antibody) and Russell's viper venom factor X activator-labeled antibody was captured onto an anti-chicken immunoglobulin G- or avidin-coated plate. The approach was found to be highly sensitive, but the assay relies on a sophisticated amplification system

utilizing a snake venom coagulation factor and is limited by its complexity and reagent expense. Poli et al. (24) developed a modified enzyme-linked immunosorbent assay (ELISA) method for serotypes E and F using affinity-purified horse polyclonal antibodies directed against the 50-kDa C fragment of both toxin E and toxin F. An extensive collaborative study for an amplified ELISA for A, B, E, and F was reported by Ferreira et al. (16). The amplified ELISA was compared with the AOAC International official method for the detection of *C. botulinum* toxins in foods and found suitable to be a preliminary test for the toxin detection and typing. Fach et al. (14) developed a PCR-based identification system which was able to detect even 1 cell of *C. botulinum* type A or B and 10 cells of *C. botulinum* type E strains in the samples. The method is sensitive; however, it can take a relatively longer time to perform (a minimum of 5 to 6 h for ELISA and overnight culture enrichments) than the biosensor for *C. botulinum* toxin detection developed by Ogert et al. (23) and Singh and Silvia (29), which offers fast throughput screening but requires expensive instrumentation. Therefore, it is important to search for cost-effective and practical assay methods capable of analyzing a large quantity of samples in the field. These monitoring assays must be robust and easy to perform.

In recent years, there has been a growing interest in developing low-cost lateral-flow assays (LFDs) for the rapid identification of analytes and pathogens (2, 19, 26, 25, 30, 31). Primarily, these assays are based on an immunochromatographic procedure that utilizes antigen-antibody reaction on a nitrocellulose membrane that is indicated by a color band from attached gold beads. The sensitivity of the lateral-flow assay (10 to 20 ng/ml) was found to be much less than those of certain ELISA kits (100 to 1,000 pg/ml). While conventional lateral-flow assays provide only a visual assay result able to

* Corresponding author. Mailing address: Center for Food Safety and Applied Nutrition, Food and Drug Administration, 5100 Paint Branch Pkwy., HFS-302, College Park, MD 20740-3835. Phone: (301) 436-1570. Fax: (301) 436-2632. E-mail: Shashi.Sharma@cfsan.fda.gov.

Report Documentation Page			Form Approved OMB No. 0704-0188		
Public reporting burden for the collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington VA 22202-4302. Respondents should be aware that notwithstanding any other provision of law, no person shall be subject to a penalty for failing to comply with a collection of information if it does not display a currently valid OMB control number.					
1. REPORT DATE 2005		2. REPORT TYPE N/A		3. DATES COVERED -	
4. TITLE AND SUBTITLE Evaluation of Lateral-Flow Clostridium botulinum Neurotoxin Detection Kits for Food Analysis				5a. CONTRACT NUMBER	
				5b. GRANT NUMBER	
				5c. PROGRAM ELEMENT NUMBER	
6. AUTHOR(S)				5d. PROJECT NUMBER	
				5e. TASK NUMBER	
				5f. WORK UNIT NUMBER	
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) U.S. Food and Drug Administration, Center for Food Safety and Applied Nutrition, College Park, Maryland 20740				8. PERFORMING ORGANIZATION REPORT NUMBER	
9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES) Naval Medical Research Center 503 Robert Grant Avenue Silver Spring, MD 20910-7500				10. SPONSOR/MONITOR'S ACRONYM(S)	
				11. SPONSOR/MONITOR'S REPORT NUMBER(S)	
12. DISTRIBUTION/AVAILABILITY STATEMENT Approved for public release, distribution unlimited					
13. SUPPLEMENTARY NOTES					
14. ABSTRACT					
15. SUBJECT TERMS					
16. SECURITY CLASSIFICATION OF:			17. LIMITATION OF ABSTRACT SAR	18. NUMBER OF PAGES 7	19a. NAME OF RESPONSIBLE PERSON
a. REPORT unclassified	b. ABSTRACT unclassified	c. THIS PAGE unclassified			

detect analytes in the nanogram range, they are user friendly, relatively inexpensive, and ideally suited for on-site testing by minimally trained personnel and can be adaptable for high-throughput laboratory or field use. However, lateral-flow assays do have a major limitation: they deliver only qualitative results, i.e., a "yes" or "no" answer detected by the human eye. Manufacturers have been developing rapid tests that deliver quantitative results, some of which are being commercialized with some success. The most common readers translate line intensity into analyte concentration and use either colorimetric reflectance or a charge-coupled-device camera to measure the signal intensity. However, like other rapid assays, sophisticated instruments limit the detection time and require trained analysts. Other commercially available readers use the emittance of a fluorescent label that is fixed to conjugate particles. In this study, we evaluated two immunodetection kits for the detection of BoNT type A (BoNT/A), BoNT/B, and BoNT/E in foods and developed procedures to extract the toxin from food samples.

MATERIALS AND METHODS

Liquid, solid, and semisolid foods were purchased from local grocery stores. Fresh raw milk was obtained from the USDA Agriculture Research Service (Beltsville, MD). Ice cream was kept at -20°C , while other food samples were kept at 4°C until used. Bovine serum albumin was purchased from Sigma Chemical Co. (St. Louis, MO). Hemagglutinin-33 was purified from *C. botulinum* type A according to the procedure of Fu et al. (17).

Neurotoxins. Pure neurotoxin complexes of types A, B, and E were purchased from Metabio, Inc. (Madison, WI). They were prepared from A strain Hall, B strain Okra, E strain Alaska, and type F strain Langeland. A sodium dodecyl sulfate-polyacrylamide gel electrophoresis gel obtained from Metabio indicated a pure neurotoxin complex preparation. The protein concentration was determined at 595 nm on an ELx801 Ultra microplate ELISA plate reader (Bio-Tek Instrument, Winooski, VT) using the Bio-Rad protein assay kit according to the procedure of Bradford (5).

***C. botulinum* detection kits.** The Bot-Tox-BTA kit for the detection of BoNT/A was purchased from Alexeter Technologies (Tetracore), and kits from BioThreat Alert (Gaithersburg, MD) and from the Naval Medical Research Center (NMRC), Silver Spring, MD, were purchased for detecting BoNT/A, -/B, and -/E. These kits can be obtained commercially. These tests utilize a combination of antibodies to selectively detect the *C. botulinum* neurotoxin in aqueous samples. The methods and protocols used to make these devices have proprietary rights reserved with the manufacturers, and the source, type of antibody, purification, characterization, and the affinity protocols have not been disclosed by the Department of Defense and Alexeter Technologies. The kits provided by the Naval Medical Research Center were developed by the U.S. Army Medical Research Institute for Infectious Diseases' Critical Reagent Program for Chemical and Biological Defense at the Aberdeen Proving Ground in MD. The kits were used according to the manufacturers' instructions.

Culture conditions and growth of *C. botulinum* strains. Type A strains were grown anaerobically in a Bacto cooked meat medium (CMM; Difco Laboratories, Detroit, MI) overnight at 35°C . One milliliter of overnight-grown culture was transferred to Trypticase-peptone-glucose-yeast extract broth (TPGY) and grown at 35°C for 5 days. The culture material was clarified by centrifugation at $10,000 \times g$ for 10 min followed by sterile filtration through a $0.45\text{-}\mu\text{m}$ -pore-size bacteriological Acrodisc filter. Culture filtrate was stored at 4°C until used.

Sample preparation. Control samples of botulinum neurotoxin complexes A, B, and E (range, 0.2 to $1\text{ }\mu\text{g/ml}$) were prepared in 500 μl of sample buffer supplied with each test kit. For liquid food samples such as orange juice, bottled water, soft drinks, vanilla extract, and apple juice, a 5-ml food sample was spiked with 100 ng/ml of pure neurotoxin complex A or B or 200 ng/ml of pure neurotoxin complex E and incubated for 30 min at room temperature (25°C). The samples were then centrifuged at $7,000 \times g$ (Eppendorf Brinkmann Inc.) for 30 min at 4°C to remove solid particles. Subsequently, 500 μl of the supernatant was thoroughly mixed with 500 μl of sample buffer in a glass test tube and used for the assay.

High-fat-content and viscous foods such as ice cream, milk, and honey were spiked with 100 ng/ml of pure botulinum neurotoxin complex A, B, or E and

incubated for 30 min at room temperature (25°C). The samples were then diluted to a ratio of 1:5 in a sample buffer (0.01 M phosphate buffer, pH 7.4). After incubation, the samples were centrifuged at $7,000 \times g$ for 30 min at 4°C to remove solid particles and/or the lipid layer. Subsequently, 500 μl of the supernatant was thoroughly mixed with 500 μl of sample buffer in a glass test tube and used for the assay.

For semisolid food samples, 10 g of food was spiked with 200 ng/ml of botulinum neurotoxin complex A or B or 400 ng/ml of neurotoxin complex E in separate vials and incubated for 30 min at room temperature (25°C). After incubation, the samples were mixed with 10 ml of sample buffer and homogenized with a bench top stomacher (Seward, Cincinnati, OH) to make a homogeneous suspension. For solid food, a 10-g portion of the food sample was spiked with 200 ng/ml of pure neurotoxin complex A, B, or E and incubated at room temperature (25°C) for 30 min. The portion was then chopped into small pieces by hand, thoroughly mixed with 10 ml of sample buffer, and homogenized with a bench top stomacher to make a homogeneous suspension. The food-buffer mixture was centrifuged at $7,000 \times g$ for 30 min at 4°C to remove solid particles. Five hundred microliters of food sample supernatant was thoroughly mixed with 500 μl of sample buffer in a glass test tube and used for the experiment.

An organic solvent extraction was performed for milk-based food samples. A 5-ml food sample was spiked with 200 ng/ml of pure neurotoxin complex A or B or 400 ng/ml of neurotoxin complex E and incubated at room temperature (25°C) for 30 min. After incubation, the sample was centrifuged at $7,000 \times g$ for 10 min to remove the solid particles. The supernatant was carefully removed and mixed with 5 ml of 50% methanol (vol/vol) in sterilized distilled water. The sample was then thoroughly mixed and allowed to settle at room temperature (25°C). The upper layer containing solvent and lipids was removed by Pasteur pipette, and the remaining lower portion of aqueous mixture was used for the experiment.

In order to reduce the matrix interference of samples which were found negative when they were undiluted, we spiked foods with a 10-fold-higher concentration of toxin and then diluted them 10-fold in a 0.01 M phosphate buffer, pH 7.4, to adjust the final concentration to 10 ng/ml for BoNT/A and -/B and 20 ng/ml for BoNT/E. To check dilution error, the protein concentrations of undiluted and diluted samples were determined at 595 nm on an ELx801 Ultra microplate ELISA plate reader (Bio-Tek Instrument, Winooski, VT) using the Bio-Rad protein assay kit according to the procedure of Bradford (5). The standard error was calculated to be 0.0166 ($n = 3$) for diluted protein concentrations in the buffer samples. Therefore, we assumed that the diluted samples contained protein concentrations of 10 ng/ml for BoNT/A or -/B or 20 ng/ml for BoNT/E.

Principle and test procedure. The assay works on an immunochromatographic principle. The BoNT-containing liquid sample rehydrates the antitoxin antibody that is labeled with gold nanoparticles. The toxin-antibody complex then migrates via capillary action along the nitrocellulose membrane placed between two plastic strips. In the sample window, the BoNTs bind with the antitoxin antibody-gold particle conjugate and forms a red line in the sample window. The intensity of the red color on the test line is directly proportional to the amount of BoNT present in the sample. The red line on the control window acts as a positive control to assure that functional, conjugated antibody has migrated throughout the system (Fig. 1A). A colored line in the control window should always appear regardless of the presence or absence of neurotoxin in the sample. For each food sample and culture filtrate, three individual tests were performed in triplicate ($n = 3$). Unless otherwise stated, the results represent all three tests, whether positive or negative. For sensitivity studies, each toxin type (BoNT/A, -/B, or -/E) was tested in triplicate with three independent experiments ($n = 3$).

Each test device was removed from a protective pouch and placed on a flat surface. The device was labeled with control or sample identifications. Samples of 150 μl were placed into the round sample port. The manufacturer recommends 15 min for recording the results. In our experiments, results were recorded visually after 15 min and 30 min of incubation at room temperature (25°C). In the absence of BoNT, no line develops in the sample window, indicating a negative result (Fig. 1B).

RESULTS

Sensitivity. A series of control samples containing BoNT/A, -/B, and -/E (0.2 to 100 ng/ml) was tested to determine the sensitivity of the test. The sensitivity tests for BoNT/A and -/B and those for BoNT/E were performed separately. No colored lines were observed in the sample windows for samples containing 0.2 to 9 ng/ml of BoNT/A or -/B or for samples con-

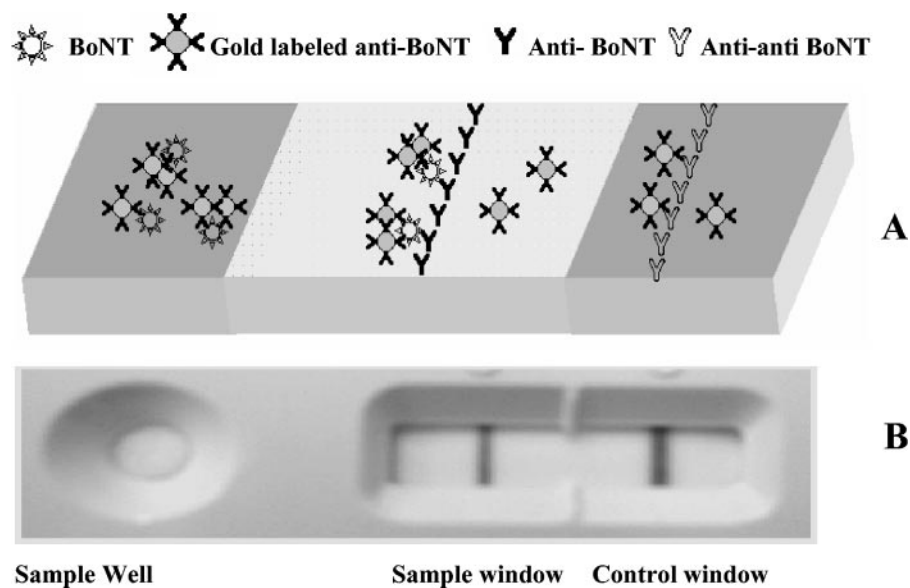


FIG. 1. (A) Schematic representation of the lateral-flow assay. A mixture of target BoNT and antibody migrates by capillary force to the membrane toward the capture zone. (B) Overview of a lateral-flow assay. The antibodies specific to a toxin immobilized onto a membrane at the test line position as shown in the sample window in the panel. These antibodies will react only with the BoNT toxin. There is also a control line, as shown in the control window in the figure, which indicates whether the test has worked or not. If a sample is positive for BoNT, there will be two lines in the viewing windows (sample and control); if negative, there will be only one at the control line position. If there is just one line at the test line position, the result should be ignored and considered invalid.

taining 0.2 to 19 ng/ml of BoNT/E. Positive red lines appeared in the sample windows for the samples containing 10 ng/ml of either BoNT/A or -/B. In the case of BoNT/E, the red line appeared with samples containing 20 ng/ml of the toxin, suggesting that the test can detect concentrations of BoNT/A and -/B as low as 10 ng/ml and of BoNT/E as low as 20 ng/ml.

Milk products. The Alexeter Technologies kit showed poor filtration and migration of the high-fat-content samples compared to that of the NMRC kit. These samples filtered poorly through the filtration device and were unable to migrate to the control window (Table 1). However, a clear red-colored line appeared in the sample window, indicating a positive reaction. With the NMRC kit, the samples, except for whipping cream, reached the sample and control windows within 15 min (Table 2). Both BoNT/A and -/B assay kits exhibited similar results. In

the case of the BoNT/E assay, samples such as whipping cream, half-and-half, and raw milk showed negative results. Very slow migrations of the samples were observed (Table 2).

To improve the filtration of the products such as pasteurized milk, 1% milk, half-and-half, and whipping cream, the samples were centrifuged at $7,000 \times g$ for 10 min. The supernatant was carefully taken out through a pipette and directly applied on the kit. In another experiment, samples were diluted to a 1:10 ratio in a sample buffer containing 0.01 M phosphate-buffered saline, pH 7.4 (final concentrations, 10 ng/ml for BoNT/A and -/B and 20 ng/ml for BoNT/E). Both centrifugation and further dilution markedly improved the filtration process, and all four samples were able to filter through the device and migrate (Tables 1 and 2).

It was assumed that the slow filtration could be due to the

TABLE 1. Detection of BoNT/A and -/B in milk samples using Alexeter Technologies-developed lateral-flow in vitro bioassay^a

Milk product ^b	BoNT/A results for samples that were:								BoNT/B results for samples that were:							
	Undiluted		Diluted 1:10		Centrifuged		Methanol extracted		Undiluted		Diluted 1:10		Centrifuged		Methanol extracted	
	15 min	30 min	15 min	30 min	15 min	30 min	15 min	30 min	15 min	30 min	15 min	30 min	15 min	30 min	15 min	30 min
1% Fat*	+	+	+	+	+	+	—	—	+	+	+	+	+	+	—	—
Fat free	+	+					—	—	+	+					—	—
Whipping cream	—	—	+	+	+	+	—	—	—	—	+	+	+	+	—	—
Milk*	+	+					—	—	+	+					—	—
Half-and-half	—	—	+	+	+	+	—	—	—	—	+	+	+	+	—	—
2% Fat	+	+					—	—	+	+					—	—
Raw milk	—	—	+	+	+	+	—	—	—	—	+	+	+	+	—	—
Control	+	+					—	—	+	+					—	—

^a The food samples were spiked with 10 ng/ml of BoNT/A or BoNT/B. The results were recorded after 15 and 30 min. The results represent the mean values of three samples for each food type ($n = 3$). No mixed results were observed.

^b An asterisk indicates that the sample reached to the test window but not to the control window.

TABLE 2. Detection of BoNT/A, -/B, and -/E in milk samples using NMRC-developed lateral-flow in vitro bioassay^a

Milk product ^b	BoNT/A results for samples that were:								BoNT/B results for samples that were:								BoNT/E results for samples that were:							
	Undiluted		Diluted 1:10		Centrifuged		Methanol extracted		Undiluted		Diluted 1:10		Centrifuged		Methanol extracted		Undiluted		Diluted 1:10		Centrifuged		Methanol extracted	
	15 min	30 min	15 min	30 min	15 min	30 min	15 min	30 min	15 min	30 min	15 min	30 min	15 min	30 min	15 min	30 min	15 min	30 min	15 min	30 min	15 min	30 min	15 min	30 min
1% Fat*	+	+					-	-	+	+					-	-	+	+	+	+	+	+	-	-
Fat free	+	+					-	-	+	+					-	-	+	+	+	+	+	+	-	-
Whipping cream	-	-	+	+	+	+	-	-	-	-	+	+	+	+	-	-	-	-	+	+	+	+	-	-
Milk*	+	+					-	-	+	+					-	-	+	+	+	+	+	+	-	-
Half-and-half	+	+	+	+	+	+	-	-	+	+	+	+	+	+	-	-	-	-	+	+	+	+	-	-
2% Fat	+	+					-	-	+	+					-	-	+	+	+	+	+	+	-	-
Raw milk	+	+					-	-	+	+					-	-	-	-	+	+	+	+	-	-
Control	+	+					-	-	+	+					-	-	+	+	+	+	+	+	-	-

^a Food samples were spiked with 10 ng/ml of BoNT/A or BoNT/B and 20 ng/ml of BoNT/E. The results were recorded after 15 and 30 min. The results represent the mean values of three samples for each food type ($n = 3$). No mixed results were observed.

^b An asterisk indicates that the sample reached to the test window but not to the control window.

presence of lipid in the samples and that an organic solvent may dissolve or separate the fat from the aqueous solution. Therefore, we used the organic solvent methanol to partition the toxin. After extraction, the lower aqueous layer was applied to the kit. No color line appeared in either the sample or the control window. Similar results were observed for both manufacturers' kits (Tables 1 and 2).

Solid foods, liquid foods, and seafood. A wide variety of food samples, including liquid, solid, and semisolid food and seafood, were tested. Of the 33 solid and liquid foods, only 10 foods showed a positive test endpoint within 15 min when used in undiluted forms. Ice cream, orange juice, and honey gave negative results with the Alexter Technologies kit. These food samples were found positive after being diluted to a ratio of 1:10 in a sample buffer (Table 3).

With the NMRC assay kit, a majority of liquid foods such as orange juice, bottled water, Coca-Cola, and vanilla extract showed positive results within 15 min. Solid foods such as broccoli, allspice, cinnamon, catfish nugget, snow crab, and Atlantic salmon also showed positive reactions within 15 min.

TABLE 3. Detection of *Clostridium botulinum* neurotoxin in food samples using Alexter Technologies' lateral-flow-based in vitro bioassay^a

Food matrix	BoNT/A				BoNT/B			
	Undiluted		Diluted		Undiluted		Diluted	
	15 min	30 min	15 min	30 min	15 min	30 min	15 min	30 min
Milk	+	+	+	+	+	+	+	+
Ice cream	-	-	+	+	-	-	+	+
Orange juice	-	-	+	+	-	-	+	+
Apple juice	+	+	+	+	+	+	+	+
Bottled water	+	+	+	+	+	+	+	+
Coca-Cola	+	+	+	+	+	+	+	+
Honey	-	-	+	+	-	-	+	+
Control	+	+	+	+	+	+	+	+

^a Only those samples which were found negative when they were undiluted were diluted to ratios of 1:10 (final concentrations after dilution in solution, 10 ng/ml for BoNT/A and -/B). The dilutions were done after spiking the samples with toxin. The results were recorded after 15 and 30 min. The results represent the mean values of three samples for each food type ($n = 3$). No mixed results were observed.

The food samples of honey and infant formula were unable to reach the sample windows and showed negative reactions in the first 15 min of the test. In the case of apple juice, although the sample was filtered through a nitrocellulose membrane, the membrane turned slightly red, making it difficult to read the result. Only those viscous samples which were found negative within 30 min in undiluted form in a sample buffer were diluted (final concentrations, 10 ng/ml for BoNT/A and -/B and 20 ng/ml for BoNT/E). All diluted food samples turned positive within 15 or 30 min of test time (Table 4).

Unlike the BoNT/A and -/B kits, in which the times to positive observation vary, the results virtually remained unchanged even after 30 min for BoNT/E. The undiluted liquid sample either did not filter or was unable to reach the control window, leaving a negative result. However, upon dilution, in a manner similar to that seen in the assays for BoNT/A and -/B, food samples were able to reach the control windows. For BoNT/E, although the test was considered positive and both sample and control window lines were clear and visible, the color was less intense than that for BoNT/A and BoNT/B (Table 4).

Detection of cultured toxin and specificity. Seventy-seven strains of *Clostridium botulinum* type A isolated from different environmental and clinical samples were chosen for the detection of the cultured toxin. These strains are isolates of out-breaks and are from a Food and Drug Administration *Clostridium botulinum* depository. The culture filtrate from each isolate was diluted 100-fold in a 0.01 M phosphate buffer, pH 7.4, before being used in the assay. Further attempts to dilute the cultured toxin (1×10^3 dilution) failed, and negative results were observed. Except for strain CS-A, all strains showed positive results when grown in TPGY medium. Culture filtrates of four strains, i.e., SKOR1-A (8833B1), SKOR@-A (8933E1), CS-A, and 8-A, grown in a CMM medium were found negative for the presence of BoNT/A. In order to test for false-positive results, we used denatured BoNT/A complex (heated at 100°C for 10 min), pure BoNT/E and -/F complexes, tap water, bovine serum albumin, purified hemagglutinin-33 (from BoNT/A complex), and culture filtrates of 26 other non-botulinum strains, including strains of *Clostridium acetobutylicum*, *Clostridium chauvoei*, *Clostridium sporogenes*, *Clostridium*

TABLE 4. Detection of *Clostridium botulinum* neurotoxin in food samples using lateral-flow-based in vitro bioassay

Food matrix	NMRC assay result for ^a :											
	BoNT/A				BoNT/B				BoNT/E			
	Undiluted		Diluted		Undiluted		Diluted		Undiluted		Diluted	
	15 min	30 min	15 min	30 min	15 min	30 min	15 min	30 min	15 min	30 min	15 min	30 min
Ice cream	—	+	+	+	—	—	+	+	—	—	+	+
Orange juice	+	+			—	+	+	+	—	—	+	+
Apple juice	—	—	+	+	—	—	+	+	+	+		
Bottled water	+	+			+	+			+	+		
Coca-Cola	+	+			+	+			+	+		
Honey	—	—	+	+	—	—	+	+	+	+		
Infant formula	—	+	+	+	—	+	+	+	—	—	+	+
Vanilla extract	+	+			+	+			+	+		
Broccoli	+	+			+	+			+	+		
Oregano	—	+	+	+	—	+	+	+	—	—	+	+
Potato salad	—	+	+	+	—	+	+	+	—	—	+	+
Allspice	+	+			+	+			—	—	+	+
Basil leaves	—	+	+	+	—	—	+	+	—	+	+	+
Cinnamon	+	+			+	+			—	+	+	+
Black pepper	—	+	+	+	—	—	+	+	—	—	+	+
Cumin	—	+	+	+	—	+	+	+	—	—	+	+
Curry powder	—	+	+	+	—	+	+	+	—	—	+	+
Ginger	—	+	+	+	—	+	+	+	—	—	+	+
Paprika	—	+	+	+	—	—	+	+	+	+		
Smoked turkey	—	+	+	+	—	+	+	+	—	—	+	+
Turkey breast	—	+	+	+	—	+	+	+	—	+	+	+
Pastrami	—	+	+	+	—	+	+	+	+	+		
Roast beef	—	+	+	+	—	+	+	+	—	+	+	+
Cured ham	—	+	+	+	—	+	+	+	+	+		
Pepperoni	—	+	+	+	—	+	+	+	—	+	+	+
Provolone	—	+	+	+	—	+	+	+	—	—	+	+
Hard salami	—	+	+	+	—	+	+	+	—	+	+	+
Cooked ham	—	+	+	+	—	+	+	+	+	+		
Domestic Swiss cheese	—	+	+	+	—	+	+	+	+	+		
Catfish nugget	+	+			+	+			+	+		
Shrimp	—	—	+	+	—	—	+	+	—	—	+	+
Snow crab leg	+	+			+	+			+	+		
Atlantic salmon	+	+			+	+			+	+		
Control	+	+			+	+			+	+		

^a Only those samples which were found negative when they were undiluted were diluted to ratios of 1:10 (final concentrations after dilutions in solution, 10 ng/ml for BoNT/A and -/B and 20 ng/ml for BoNT/E). The dilutions were done after spiking the samples with toxin. The results were recorded after 15 and 30 min. The results represent the mean values of three samples for each food type ($n = 3$). No mixed results were observed.

histolyticum, *Clostridium tetani*, *Clostridium bifermentans*, *Clostridium sordellii*, *Clostridium perfringens*, *Bacillus subtilis*, *Bacillus cereus*, and *Bacillus megaterium* as interfering agents. No false positives were observed for 26 of the culture filtrates of 26 strains or for other interfering agents except for BoNT/F. A very faint positive line appeared in sample windows with BoNT/F samples.

DISCUSSION

Immunoassay technologies are ideal for the qualitative and quantitative detection of many types of proteins and pathogens in complex matrices (1, 6). Lateral-flow assays have been used extensively as diagnostic tools for monitoring toxins (7, 32, 33). However, their effectiveness is confirmed only for toxins in serum and environmental samples. We wanted to evaluate lateral-flow assays as an early warning tool for the detection of *C. botulinum* toxins in a variety of foods. In this investigation, two assay kits appear to be similar in technology but performed slightly differently in our experiments. The major difference

observed seems to be in the mechanical assemblies of the filtration devices and the nitrocellulose membrane supports. For BoNT/A, the assay sensitivities were found to be similar for both the Alexeter Technologies assay kit and the NMRC assay kit. During the cross-reaction studies, we observed a positive reaction in the BoNT/A kit with BoNT/B neurotoxin. When tested for sensitivity, the BoNT/A kit was found to be equally sensitive for BoNT/B neurotoxin. The cross-reaction between the BoNT/A and BoNT/B antibodies has been reported in the literature (8, 10). Alexeter Technologies, Inc. markets kits only for BoNT/A. The sensitivity limit of the NMRC BoNT/E kit was found to be 20 ng/ml. The assay's sensitivities for BoNT/A, -/B, and -/E were much less than those of the mouse bioassay and therefore, for the comparison, we did not perform the mouse bioassay.

The assays were visualized as red lines created by the bound colloidal gold; therefore, sensitivity is limited to what can be seen by the human eye and is lower than that of sensitive instrumentation such as a spectrophotometer or a fluorescent reader. Typically, an arbitrary quantitation of the detection

TABLE 5. Milk product contents^a

Milk product ^b	Serving size (ml)	Total fat (g)	Saturated fats (g)	Cholesterol (mg)	Sodium (mg)	Carbohydrates (g)	Dietary fiber (g)	Sugars (g)	Proteins (g)
1% Low fat	50	0.52	0.31	2.08	27.08	2.50	0.00	2.29	1.67
Fat free	50	0.00	0.00	1.04	27.08	2.50	0.00	2.29	1.67
Whipping cream	50	16.67	11.67	66.67	16.67	0.00	0.00	0.00	0.00
Vitamin D	50	1.67	1.04	7.29	26.04	2.50	0.00	2.29	1.67
Half and half	50	5.00	3.33	25.00	25.00	3.33	0.00	1.67	1.67
2% Reduced fat	50	1.04	0.63	4.17	26.04	2.50	0.00	2.29	1.67
Raw milk	50	1.85	0.935	5.00	25.00	2.20	0.00	2.63	1.75

^a Values are normalized to 50-ml serving sizes for uniform representation.

^b Product descriptions are based on the labels attached to the products.

sensitivities of these assays is done by assigning “plus” and “minus” values, with the increasing intensity of the red line assigned a plus value. Besides the somewhat qualitative nature of this process, the weak-positive-result values can vary based on the skill of the technician responsible for validating a given lot of assays.

The lethal dose of botulinum toxin for humans is not precisely known. The estimated lethal amount of crystalline type A toxin for a 70-kg human would be approximately 0.09 to 0.15 µg intravenously or intramuscularly, 0.70 to 0.90 µg inhalationally, or 70 µg orally (4). The detection limit of the mouse bioassay is 10 pg/ml of pure BoNT/A (11). We found that the BoNT/A and -B assay kits can readily detect 10 ng/ml, while the BoNT/E kit can detect as low as 20 ng/ml of BoNT/E. Although these assays exhibited sensitivities less than the detection limit of the mouse bioassay, they are capable of detecting the toxin concentration that can cause botulism disease symptoms in humans. Therefore, they can be considered for large-scale or presumptive test screening to be followed by a confirmatory mouse bioassay.

Certain milk-based samples, such as whipping cream, half-and-half, and raw milk, when tested in undiluted form, exhibited low filtration or sometimes virtually no filtration with these test devices. Such problems were particularly encountered by the Alexeter Technologies kits. The test kits from NMRC performed relatively well in terms of filtration and sample migration of the analytes. They developed filtration and sample migration problems only with whipping cream. This was perhaps due to the high fat contents of these milk products (Table 5). In general, the filtration system prevents the solid particles and passes analytes (toxin) freely to the membrane. Lateral-flow assays work on the basis of a capillary process (2, 25) and, because of the action of capillary forces, the toxin has to migrate up the membrane to where ligands are immobilized in the capture zones. It seems that the high fat contents of the samples prevent the capillary flow of the samples. However, several other factors could cause negative results, including the inhibition of antigen antibody reaction due the presence of inhibitory enzymes present in the samples.

The manufacturers' recommended time for recording the result is 15 min. We considered an extra 15 min of incubation, because undiluted high-fat-content foods were unable to reach the control window within 15 min. In almost all cases, an extra 15 min was sufficient for the toxins to reach the control windows even with slow migration of the samples. Alternatively, dilution or centrifugation may also speed up the migration

process observed in certain milk-based samples (Table 4). Both sample and control lines appeared, indicating a positive reaction. Therefore, we recommend that at least a 30-min test time be used for high-fat-content food samples. The Alexeter Technologies kit instructions indicate that colored lines that appear after 20 min are not valid and should be ignored. We also observed slow filtration for some semisolid and solid food samples similar to that of the milk-based samples. Since these foods do not contain amounts of fat sufficient to restrict the flow of the samples through membrane and their flows were virtually similar to those of milk-based samples, it is apparent that fat is not the only factor that is responsible for slow filtration.

These lateral-flow tests can detect up to 100-fold dilutions of the cultured toxin in both TPGY medium and CMM, indicating that undiluted cultured toxins of foods can be tested for a qualitative identification of BoNT. The choice of medium and its constituents can also affect the overall results of the assay. For example, of the 77 strains grown in CMM, the assay showed three false-negative results. In the TPGY medium, no false negatives were observed, and all strains (except CS-A) showed positive reactions. Except for BoNT/F, other potentially interfering agents did not give false-positive results, perhaps indicating that a weak cross-reaction of BoNT/F antibodies occurs with BoNT/A and -B.

The ELISA method, which can detect from 0.1 to 1 ng/ml (~10 to 100 MLD) of the toxin in a variety of samples, including food, has been widely used by researchers (15) and is commercially available from MetabioLogics, Inc. The modified ELISA was compared with the AOAC International method and proved useful for screening a large quantity of samples. The assay takes approximately 5 h to complete. The major disadvantage of ELISA is that it is relatively time and labor intensive. Unlike the typical ELISA, in which either antibodies or antigens are directly bound to a polystyrene or polyvinyl microplate, electrochemiluminescence detection offers relatively better sensitivity (0.05 to 1.0 ng/ml, ~10 MLD) through a process in which antibodies are attached on the magnetic particles in a large volume of sample suspension and then captured and bound in a small area. Electrochemiluminescence detection provides increased sensitivity due to high luminescent-signal-to-noise ratios. However, the assay performance could vary significantly depending on the sample matrices. Because of this, matrix-specific positive and negative control samples are used to establish standard curves and cutoff values. An expensive instrumentation is also a major limitation with

the assay. Lateral-flow assays are simple to use and require minimal training. They are relatively cost effective and can be used as first-responder detection kits. Although these devices sound ideal, they have limitations. Key limitations are that they are less sensitive than other in vitro detection systems and that the assessments of results are strictly qualitative. New-generation lateral-flow assays based on fluorescence have substantially greater sensitivities and dynamic ranges than did the ones used in this study. Nevertheless, the development of such assays has been limited by the need for stable dyes that do not cause sample interference and by the fact that the instrumentation required was both complex and expensive.

In summary, the performance of the LFDs indicates that these assays are suitable as screening tools for a variety of foods potentially contaminated with BoNT. The speed of the assay (15 min), in concert with proven precision, demonstrates its usefulness in testing samples when used by trained laboratory operators. However, an extra 15 min of incubation time needs careful evaluation and enzyme and antigen antibody standardization before it can be used for large-scale screening. However, due to their relatively insensitive natures, they can be used as presumptive tests only, and the mouse bioassay should be used for confirmatory results. Ease of operation and consistent results suggest that the LFDs can be effectively used by nontechnical staff for screening food samples.

REFERENCES

- Ahmed, F. E. 2002. Detection of genetically modified organisms in foods. *Trends Biotechnol.* **20**:215–223.
- Al-Yousif, Y., J. Anderson, C. Chard-Bergstrom, and S. Kapil. 2002. Development, evaluation, and application of lateral-flow immunoassay (immunochromatography) for detection of rotavirus in bovine fecal samples. *Clin. Diagn. Lab. Immunol.* **9**:723–724.
- Anne, C., F. Cornille, C. Lenoir, and B. P. Roques. 2001. High-throughput fluorogenic assay for determination of botulinum type B neurotoxin protease activity. *Anal. Biochem.* **291**:253–261.
- Arnon, S. S., R. Schechter, T. V. Inglesby, D. A. Henderson, J. G. Bartlett, M. S. Ascher, E. Eitzen, A. D. Fine, J. Hauer, M. Layton, S. Lillibridge, M. T. Osterholm, T. O'Toole, G. Parker, T. M. Perl, P. K. Russell, D. L. Swerdlow, and K. Tonat. 2001. Botulinum toxin as a biological weapon: medical and public health management. *JAMA* **285**:1059–1070.
- Bradford, M. M. 1976. A rapid and sensitive method for the quantitation of microgram quantities of protein utilizing the principle of protein-dye binding. *Anal. Biochem.* **72**:248–254.
- Brett, G. M., S. J. Chambers, L. Huang, and M. R. A. Morgan. 1999. Design and development of immunoassays for detection of proteins. *Food Control* **10**:401–406.
- Cangel, P., Y. Chan, K. W. Sum, K. Y. Cheung, J. F. C. Glatz, J. E. Sanderson, A. Hempel, M. Lehmann, I. Renneberg, and R. Renneberg. 2003. Development of a quantitative lateral-flow assay for rapid detection of fatty acid-binding protein. *J. Immunol. Methods* **279**:91–100.
- Chiao, D. J., R. H. Shyu, C. S. Hu, H. Y. Chiang, and S. S. Tang. 2004. Colloidal gold-based immunochromatographic assay for detection of botulinum neurotoxin type B. *J. Chromatogr. B* **809**:37–41.
- Dezfulian, M., C. L. Hatheway, R. H. Yolken, and J. G. Bartlett. 1984. Enzyme-linked immunosorbent assay for detection of *Clostridium botulinum* type A and type B toxins in stool samples of infants with botulism. *J. Clin. Microbiol.* **20**:379–383.
- Dezfulian, M., and J. G. Bartlett. 1984. Detection of *Clostridium botulinum* type A toxin by enzyme-linked immunosorbent assay with antibodies produced in immunologically tolerant animals. *J. Clin. Microbiol.* **19**:645–648.
- Doellgast, G. J., M. X. Triscott, G. A. Beard, and J. D. Bottoms. 1994. Enzyme-linked immunosorbent assay-enzyme-linked coagulation assay for detection of antibodies to *Clostridium botulinum* neurotoxins A, B, and E and solution-phase complexes. *J. Clin. Microbiol.* **32**:851–853.
- Doellgast, G. J., J. E. Brown, J. A. Koufman, and C. L. Hatheway. 1997. Sensitive assay for measurement of antibodies to *Clostridium botulinum* neurotoxins A, B, and E: use of hapten-labeled-antibody elution to isolate specific complexes. *J. Clin. Microbiol.* **35**:578–583.
- Doellgast, G. J., M. X. Triscott, G. A. Beard, J. D. Bottoms, T. Cheng, B. H. Roh, M. G. Roman, P. A. Hall, and J. E. Brown. 1993. Sensitive enzyme-linked immunosorbent assay for the detection of *Clostridium botulinum* neurotoxins A, B, and E using signal amplification via enzyme-linked coagulation assay. *J. Clin. Microbiol.* **31**:2402–2409.
- Fach, P., M. Gibert, R. Griffais, J. P. Guillou, and M. R. Popoff. 1995. PCR and gene probe identification of botulinum neurotoxin A-, B-, E-, F-, and G-producing *Clostridium* spp. and evaluation in food samples. *Appl. Environ. Microbiol.* **61**:389–392.
- Ferreira, J. L., S. Eliasberg, M. A. Harrison, and P. Edmonds. 2001. Detection of preformed type A botulinum toxin in hash brown potatoes using the amplified ELISA. *J. Assoc. Off. Agric. Chem. Int.* **84**:1460–1464.
- Ferreira, J. L., S. Maslanka, E. Johnson, and M. Goodnough. 2003. Detection of botulinum neurotoxins A, B, E, and F by amplified enzyme-linked immunosorbent assay: collaborative study. *J. Assoc. Off. Agric. Chem. Int.* **86**:314–331.
- Fu, F. N., S. K. Sharma, and B. R. Singh. 1998. A protease-resistant novel hemagglutinin purified from type A *Clostridium botulinum*. *J. Protein Chem.* **17**:53–60.
- Johnson, E. A. 1999. Clostridial toxins as therapeutic agents: benefits of nature's most toxic proteins. *Annu. Rev. Microbiol.* **53**:551–575.
- Ketema, F., C. Zeh, D. C. Edelman, R. Saville, and N. T. Constantine. 2001. Assessment of the performance of a rapid, lateral flow assay for the detection of antibodies to HIV. *J. Acquir. Immune Defic. Syndr.* **27**:63–70.
- McGrath, S., J. S. Dooley, and R. W. Haylock. 2000. Quantification of *Clostridium botulinum* toxin gene expression by competitive reverse transcription-PCR. *Appl. Environ. Microbiol.* **66**:1423–1428.
- Middlebrook, J. L. 1989. Cell surface receptors for protein toxins, p. 95–119. In L. L. Simpson (ed.), *Botulinum neurotoxin and tetanus toxin*. Academic Press, New York, N.Y.
- Niemann, H. 1991. Molecular biology of the clostridial neurotoxins, p. 303–348. In J. E. Alouf and J. H. Freer (ed.), *A sourcebook on bacterial protein toxins*. Academic Press, London, United Kingdom.
- Ogert, R. A., J. E. Brown, B. R. Singh, L. C. Shriver-Lake, and F. S. Ligler. 1992. Detection of *Clostridium botulinum* toxin A using a fiber optic biosensor. *Anal. Biochem.* **205**:306–312.
- Polli, M. A., V. R. Rivera, and D. Neal. 2002. Development of sensitive colorimetric capture ELISAs for *Clostridium botulinum* neurotoxin serotypes E and F. *Toxicon* **40**:797–802.
- Qian, S., and H. H. Bau. 2003. A mathematical model of lateral flow bioreactions applied to sandwich assays. *Anal. Biochem.* **322**:89–98.
- Rong-Hwa, S., S. Huey-Fen, L. Hwan-Wun, and T. Shiao-Shek. 2002. Colloidal gold-based immunochromatographic assay for detection of ricin. *Toxicon* **40**:255–258.
- Shone, C. C., P. Wilton-Smith, N. Appleton, P. Hambleton, N. K. Modi, S. Gatley, and J. Melling. 1985. Monoclonal antibody-based immunoassay for type A *Clostridium botulinum* toxin is comparable to the mouse bioassay. *Appl. Environ. Microbiol.* **50**:63–67.
- Singh, B. R. 2000. Intimate details of the most poisonous poison. *Nat. Struct. Biol.* **7**:617–619.
- Singh, B. R., and M. A. Silvia. 1996. Detection of botulinum neurotoxins using optical fiber-based biosensor, p. 499–508. In B. R. Singh and A. Tu (ed.), *Natural toxins II: structure, mechanism of action and detection*. Plenum Press, New York, N.Y.
- Smits, H. L., C. K. Eapen, S. Sugathan, M. Kuriakose, M. H. Gasem, C. Yersin, D. Sasaki, B. Pujianto, M. Vestering, T. H. Abdoel, and G. C. Gussenhoven. 2001. Lateral-flow assay for rapid serodiagnosis of human leptospirosis. *Clin. Diagn. Lab. Immunol.* **8**:166–169.
- Wittmann, C., and P. Y. Schreiter. 1999. Analysis of terbutylazine in soil samples by two test strip immunoassay formats using reflectance and luminescence detection. *J. Agric. Food Chem.* **47**:2733–2737.
- Wittmann, C., U. Bilitewski, T. Giersch, U. Kettling, and R. D. Schmid. 1996. Development and evaluation of a dipstick immunoassay format for the determination of atrazine residues on-site. *Analyst* **121**:863–869.
- Yeoh, H. H., and F. Sun. 2001. Assessing cyanogen content in cassava-based food using the enzyme-dipstick method. *Food Chem. Toxicol.* **39**:649–653.